

SMART PLASMA THAWER

WPFD2/4 & WPFD3/6

FAST, SAFE, ADVANCED THAWING

Rapid Dry Bath Plasma Thawer from KW Apparecchi offers safe, uniform, complete thawing of Fresh Frozen Plasma (FFP) in 15 minutes.

Available in two pocket (up to 4 plasma packs) or three pocket versions (up to 6 plasma packs)

Futureproof technology offers ease of use, total traceability and reporting of user, product and thawing via an Infra Red Barcode reader.

USB/ WiFi/ Ethernet connectivity allows the data to be easily exported.

Completely independent pockets can be used as separate thawers, starting a second thaw while already thawing in another pocket.

SPEED

- Complete thawing in 15 minutes

EASE OF USE

- Touch screen control

CONNECTIVITY

- USB/ WiFi/ SD Card/ Ethernet

TECHNICALLY ADVANCED

- Barcode reader, futureproof monitoring and reporting

FLEXIBILITY

- Independent pockets and thawing

The inbuilt thawing tracer software is fully downloadable and shows data such as bag ID number, start and end time of thaw, temperature graph and user ID.

This data is easily exported and can either be stored electronically or printed and added to patient notes providing total traceability.

Contact us to discuss your requirements and specification or to request further information.

Office


info@qedscientific.com

01663 735494

THAWING TRACER

The easy to use thawing tracer software exports all the data of a cycle and displays as a graph detailing;

- Bag Code
- Product Code
- User ID
- Start/ Finish Time
- Cycle Graph

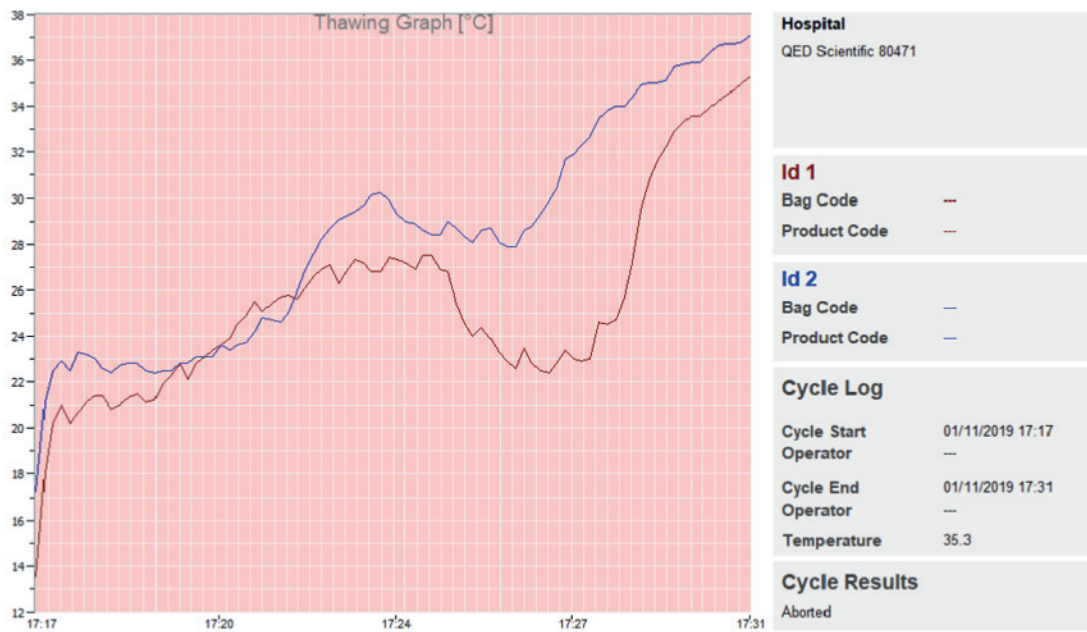

Hospital	
QED Scientific 80471	
Id 1	
Bag Code	--
Product Code	--
Id 2	
Bag Code	--
Product Code	--
Cycle Log	
Cycle Start	01/11/2019 17:17
Operator	--
Cycle End	01/11/2019 17:31
Operator	--
Temperature	35.3
Cycle Results	
Aborted	

TECHNICAL SPECIFICATION

Bag Capacity	2 x 1000ml or 4 x 460ml
User Parameters	Cycle Time Temperature Bag Mode
Screen Type	7" LCD TFT Colour Touchscreen
Service Parameters Safety	Password protected
Alarms (Audible/ Visual)	Open Door Cycle End Liquid Presence
Alarm Silence	From Touchscreen
Data Download	WiFi, USB or Ethernet

PHYSICAL SPECIFICATION

Thawer Dimensions (WxDxH)	396 x 640 x 420mm
Thawer Weight	20kg
Water Capacity	25l

STANDARDS

Classifications	Class I
Quality Standards	ISO13485, 93/42/EEC and CE Certificate
Warranty	1 year

